

The iQ Labs Floor-Mounted Series walk-in fume hood is offered in a variety of widths, depths and heights and is designed to be easily customized to fit the customers specific application. With the ability to add pass-thru openings in the side walls, stainless steel access doors, and either angled or flat front posts, the Floor-Mounted Series provides you with the versatility, durability and quality you also find in our benchtop models. With the easy to remove front panel, these hoods are perfect for chemical processes that require the need to move large pieces of a equipment in and out of the containment area.




PERFORMANCE AND AFFORDABILITY - WITHOUT COMPROMISE

- 16g cold-rolled steel superstructure
- Stainless steel sash frame and guides
- Angled posts for smooth airflow transition
- Three-piece baffle system with dampers
- LED (or explosion-proof) lighting for improved visibility
- Chemical resistant poly resin liner
- Chemical resistant powder coat finish
- ASHRAE-110 tested and approved
- Fully customizable width, height and depth option

OPTIONAL EQUIPMENT


- Remote -controlled plumbing fixtures with chrome handles, pre-piped with flex hose or hard copper
- Multiple sash options - vertical rising (standard), horizontal, or combination
- Custom pass-thru openings upon request
- Sash Stops
- Air Monitor/Alarm
- Ceiling Enclosure
- Stainless steel ramp for easy install or removal of equipment
- Blower

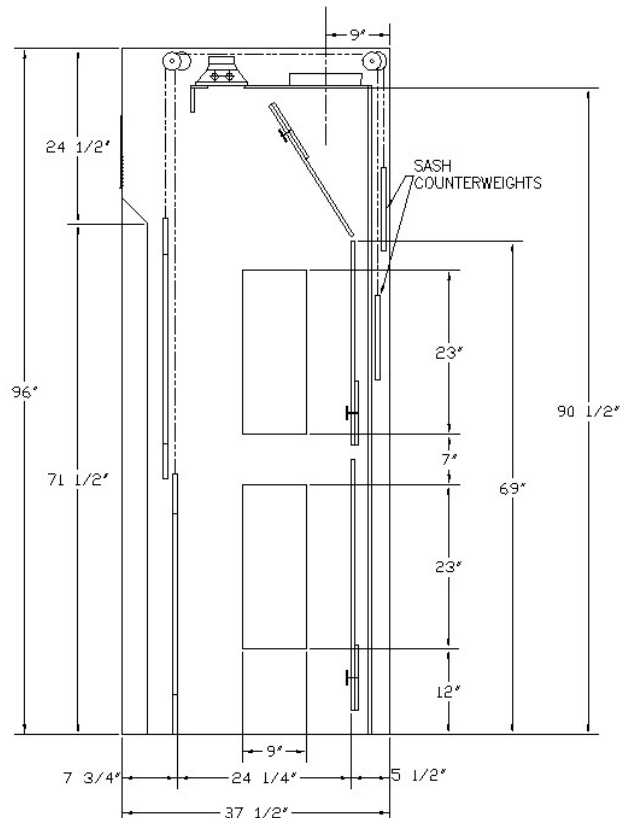

Below are the standard floor-mounted hood sizes. Custom sizes to fit specific applications are available. Please contact us directly to determine the best fit for your project and to get specific dimensioning and CFM specifications.

PART NUMBER	DIMENSIONS			INT. WORKING AREA			HOOD OPENING		CFM @ 18"*	CFM @ 32"*	HOOD WEIGHT
	W	H	D	W	H	D	W	H			
LAWI4896	48"	96"	37.5"	39.5"	69"	24.25"	39.5"	72"	493	877	925
LAWI6096	60"	96"	37.5"	51.5"	69"	24.25"	51.5"	72"	643	1144	1125
LAWI7296	72"	96"	37.5"	63.5"	69"	24.25"	63.5"	72"	793	1411	1325
LAWI8496	84"	96"	37.5"	75.5"	69"	24.25"	75.5"	72"	943	1633	1725
LAWI9696	96"	96"	37.5"	87.5"	69"	24.25"	87.5"	72"	1093	1944	1930

*CFM calculated using 100 fpm

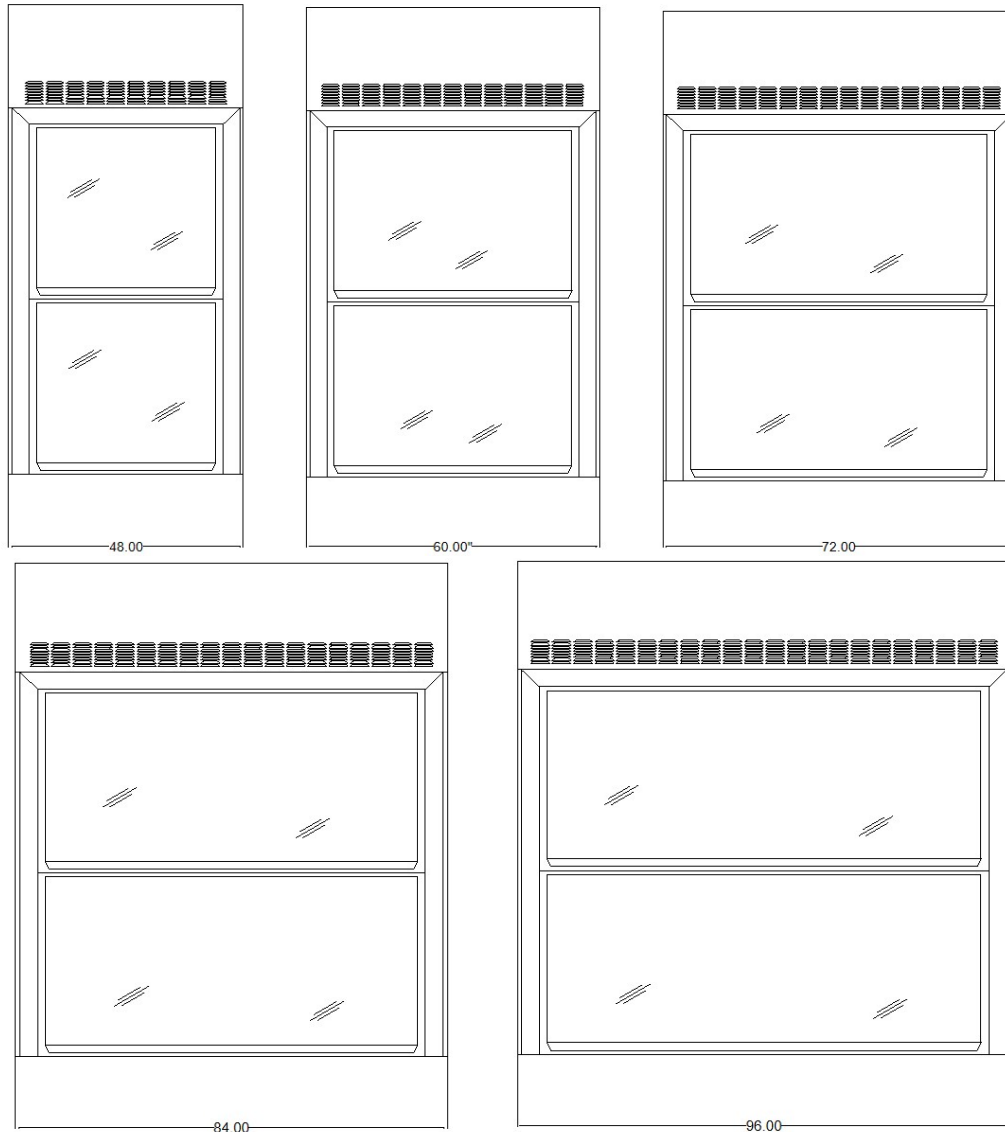

STANDARD COLOR OPTIONS (CUSTOM COLOR-MATCH AVAILABLE)

